

## SPECIFICATIONS

### Bio-Board™ Disposable Dissection Board


The **Bio-Board™ Disposable Dissection Board**, is an economical cur resistant disposable foam core board with light blue color work surface to provide color contrast against specimens. Bio-Boards™ are intended for use in routine grossing, frozen sections and photography in a professional or laboratory setting. The Disposable Dissection Board is made from Expanded polyethylene or EPE. EPE is non-crosslinked, closed-cell polyethylene foam made from low-density polyethylene resins. The board is reinforced with paper to provide a rigid workspace and will not dull sharp instruments.

- Harmonized Tariff Schedule Code: 3921.90.4010
- Country of Origin: Taiwan

### General direction of usage:

② The Bio-Board™ Disposable Dissection Boards are for single use through trained personnel only. No account should re-use Bio-Boards™ and after usage, they should be disposed as potential infective or biohazardous waste according to local laws and ordinances.

⊗ The Bio-Board™ Disposable Dissection Boards may not be used when damage of the package is shown for signs posterior breakage of the EPE foam. This creates an uneven work surface. Any damaged Bio-Boards™ should be disposed of.

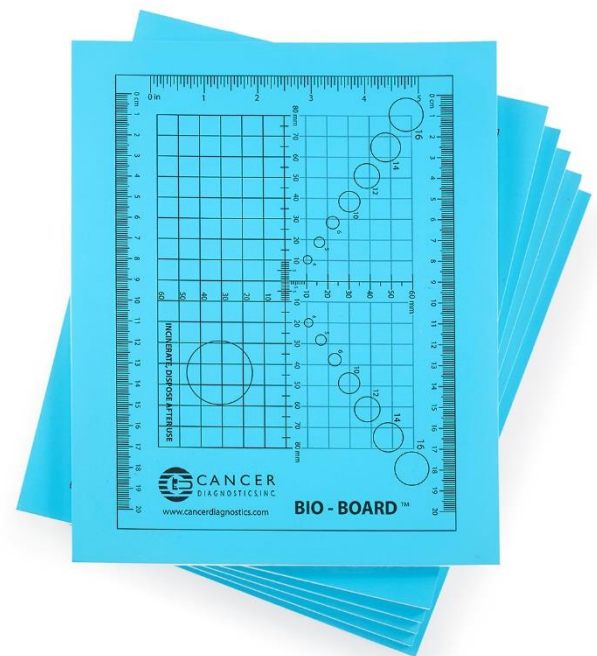

### Dimensions (in inches)

Worksurface ■		Overall ■	
<b>Length</b>	8.62	<b>Length</b>	9.875
<b>Width</b>	6.25	<b>Width</b>	7.9375
		<b>Height (Thickness)</b>	0.213

